

Thermo Scientific picoSpin 45 NMR Spectrometer

Accessible and affordable NMR spectroscopy

The Thermo Scientific™ picoSpin™ 45 spectrometer provides the power of Nuclear Magnetic Resonance (NMR) spectroscopy in a compact and affordable instrument. By drastically lowering cost and size, the picoSpin 45 instrument makes NMR spectroscopy accessible to a variety of laboratories.



The picoSpin 45 spectrometer is simple to operate. It is easy for students and technicians with limited NMR experience to use the power of NMR spectroscopy to identify chemical compounds or verify their structure. The instrument's fluid capillary system is contained within a replaceable cartridge and only requires 30 microliters of a liquid sample.

Its temperature-controlled permanent magnet does not require cryogenics, eliminating the need for consumables or custom laboratory facilities. Additionally, because the instrument is lightweight (less than five kilograms), the picoSpin 45 spectrometer is easily shared across multiple laboratories.



Thermo
SCIENTIFIC

Product Specifications

Performance

- Larmor Frequency
45 ± 1 MHz Proton
- Chemical Shift Resolution
40 ppb FWHM or Better

Magnet

- Rare Earth Permanent
- 1.1 Tesla
- Temperature Stabilized

Capillary Cartridge

- Capillary Sample Loop
400 µm Capillary ID
30 µL Sample Volume
- 1/16" Interface Fittings
- User Replaceable



Sample

- Sample Type
Liquid
- Sample Introduction
Inject Sample via Syringe

Environmental

- No Cryogenics or Compressed Gases
- Standard Laboratory Environment

Interface

- Ethernet
Direct Connection or LAN
- Embedded Web Server

Software/Collect





- Web Browser User Interface
- Data Export
JCAMP-DX
Readily opened by many NMR programs
- Data Analysis Software
Mestrelab Mnova™ NMR Software
One Year License

Dimensions

- Size: 20.3 cm × 14.6 cm × 29.2 cm
(8" × 5.75" × 11.5")
- Weight: 4.8 kg (10.5 lbs)

Electrical Requirements

- 100-230 VAC, 50/60 Hz, 150 W

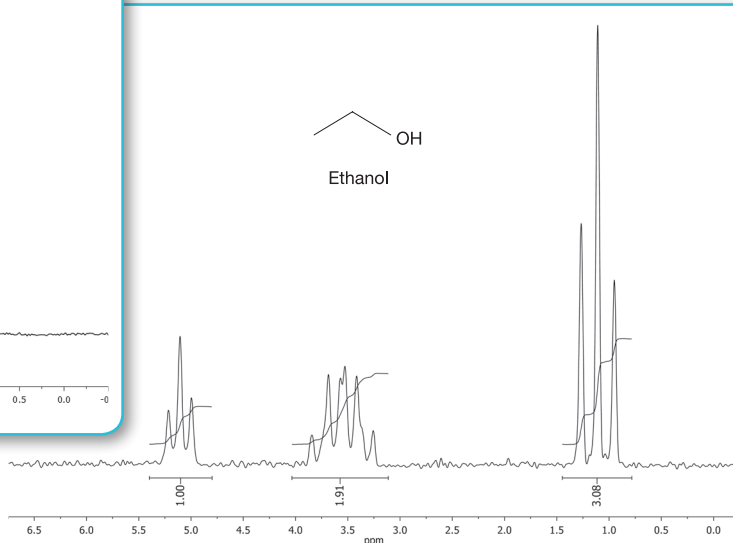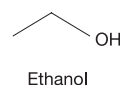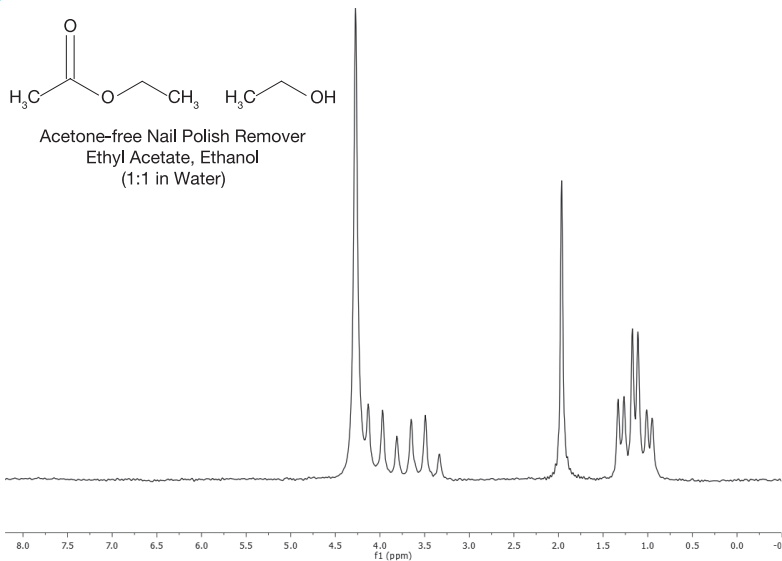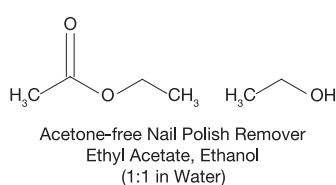

www.thermoscientific.com

©2013 Thermo Fisher Scientific Inc. All rights reserved. Mnova is a trademark of Mestrelab Research, Inc. All other trademarks are the property of Thermo Fisher Scientific and its subsidiaries. Specifications, terms and pricing are subject to change. Not all products are available in all countries. Please consult your local sales representative for details.

Africa +27 11 822 4120
Australia +61 3 9757 4300
Austria +43 1 333 50 34 0
Belgium +32 53 73 42 41
Canada +1 800 530 8447
China +86 10 8419 3588

Denmark +45 70 23 62 60
Europe-Other +43 1 333 50 34 0
Finland/Norway/Sweden
 +46 8 556 468 00
France +33 1 60 92 48 00
Germany +49 6103 408 1014

India +91 22 6742 9434
Italy +39 02 950 591
Japan +81 45 453 9100
Latin America +1 561 688 8700
Middle East +43 1 333 50 34 0
Netherlands +31 76 579 55 55

New Zealand +64 9 980 6700
Russia/CIS +43 1 333 50 34 0
Spain +34 914 845 965
Switzerland +41 61 716 77 00
UK +44 1442 233555
USA +1 800 532 4752

Thermo
 SCIENTIFIC

Part of Thermo Fisher Scientific